



**Epix<sup>®</sup>**

Electrosurgical Probes  
with Smoke Evacuation

# Epix

## Electrosurgical Probes with Smoke Evacuation

### Probe Tips

Provide precise delivery of monopolar energy



CW001



CW002



CW003




### Angled Shaft

Improves instrument reach and visibility

Reduces instrument clashing

### Built-In Smoke Evacuation

Evacuates smoke at the push of a button while minimizing loss of insufflation

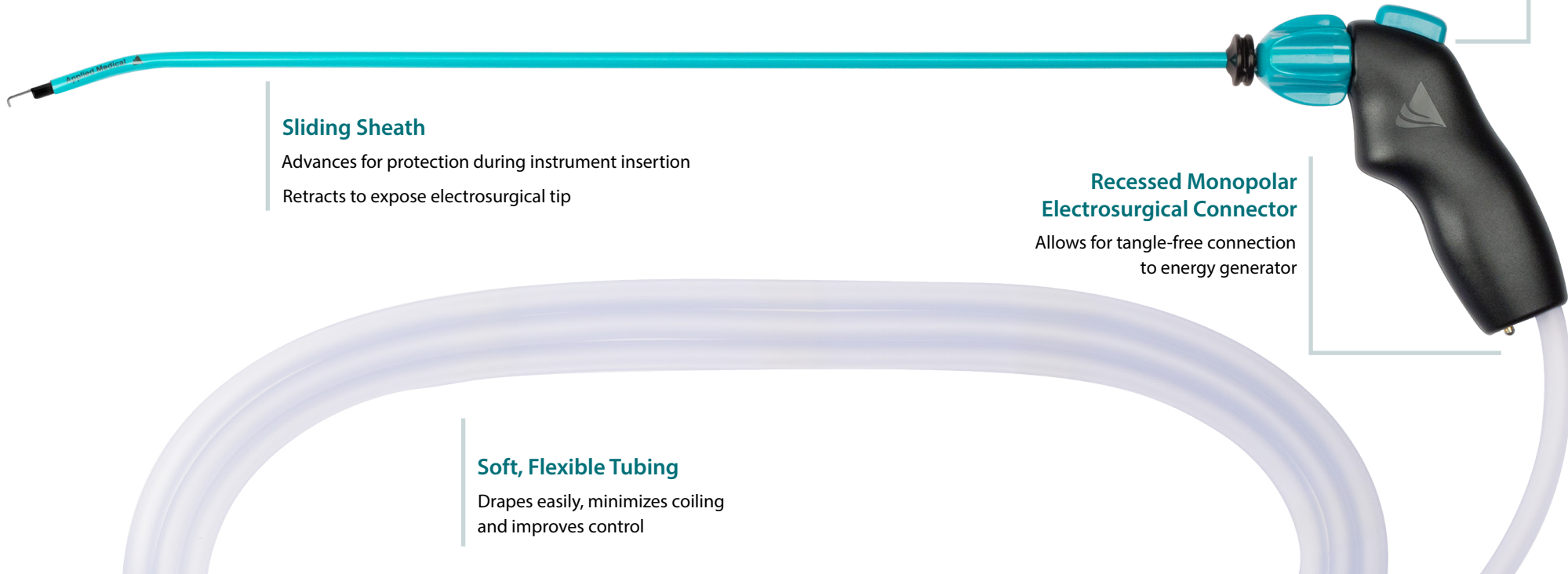

### Sliding Sheath

Advances for protection during instrument insertion

Retracts to expose electrosurgical tip

### Recessed Monopolar Electro-surgical Connector

Allows for tangle-free connection to energy generator

### Soft, Flexible Tubing

Drapes easily, minimizes coiling and improves control



## Epix Electrosurgical Probes with Smoke Evacuation

Reorder No.	Tip Design	Size	Qty/Box
CW001*	Angled Spatula	5mm x 42cm	1/Box
CW002*	Angled L-Hook	5mm x 42cm	1/Box
CW003*	Angled Needle	5mm x 42cm	1/Box

\* Models CW001, CW002 and CW003 are compatible with Applied Medical Kii 12x60-100mm and 15x100-150mm trocars.

Visit [www.appliedmedical.com/epix](http://www.appliedmedical.com/epix)  
for more information

Devices listed may not be approved in all markets.

Please contact your field implementation team member for more information on availability.

© 2019 Applied Medical Resources Corporation. All rights reserved.

Applied Medical, the Applied Medical logo design and marks designated with a ® are trademarks of Applied Medical Resources Corporation, registered in one or more of the following countries: Australia, Canada, Japan, South Korea, the United States, and/or the European Union. All of the other trademarks used herein are the property of their respective owners. 159476-EN-USA-B

